

NARROW IMPLANT

BONE LEVEL Ø2.9 mm



implantswiss

IMPLANT TECHNOLOGY

SMART
CHOICE




NARROWIMPLANT

It is developed for use in the anterior region, where the mesio-distal distance is not sufficient, or in extremely thin crest. It can be preferred in aesthetically challenging situations.

THE NARROW IMPLANT Ø 2.9 mm IS ONLY INDICATED FOR:

Single tooth restorations:

- Central & lateral incisors in the mandible
- Lateral incisors in the maxilla

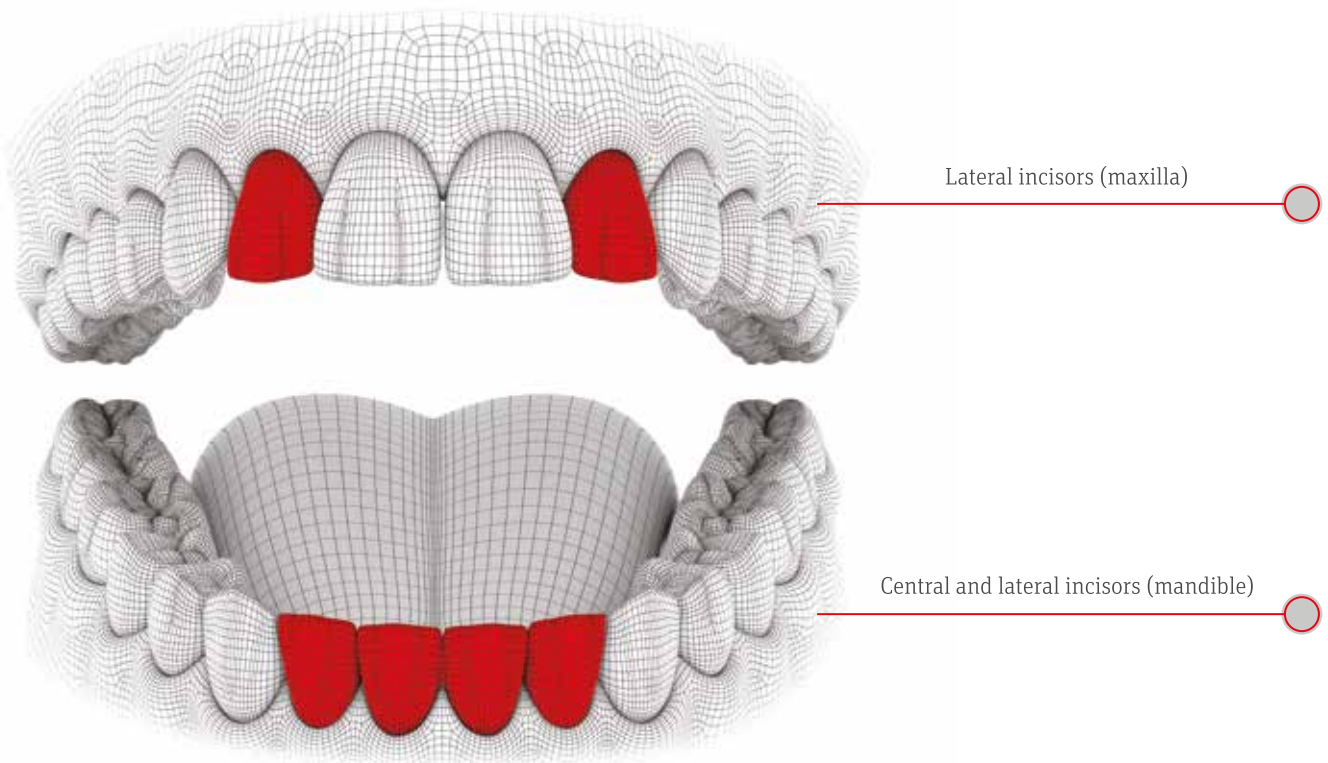


CONTRAINDICATED FOR:

- Bridges
- Edentulous


SPECIFICATIONS

- Hybrid body design
- Narrow inner connection
- No stress on the crestal bone, due to the passive thread neck design
- The threads deepening towards its apical provide compression from the apex

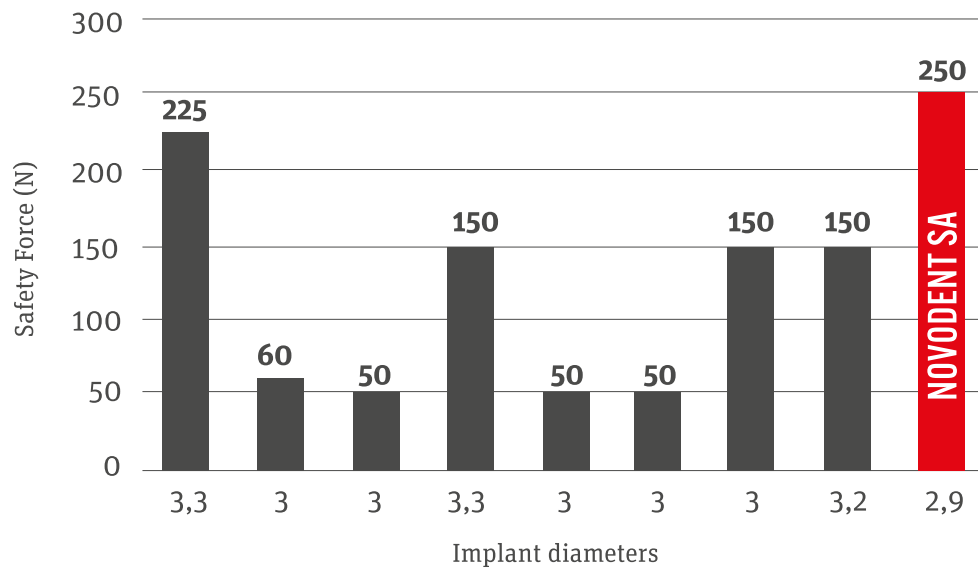

Maximum loading values for similar product groups of different brands



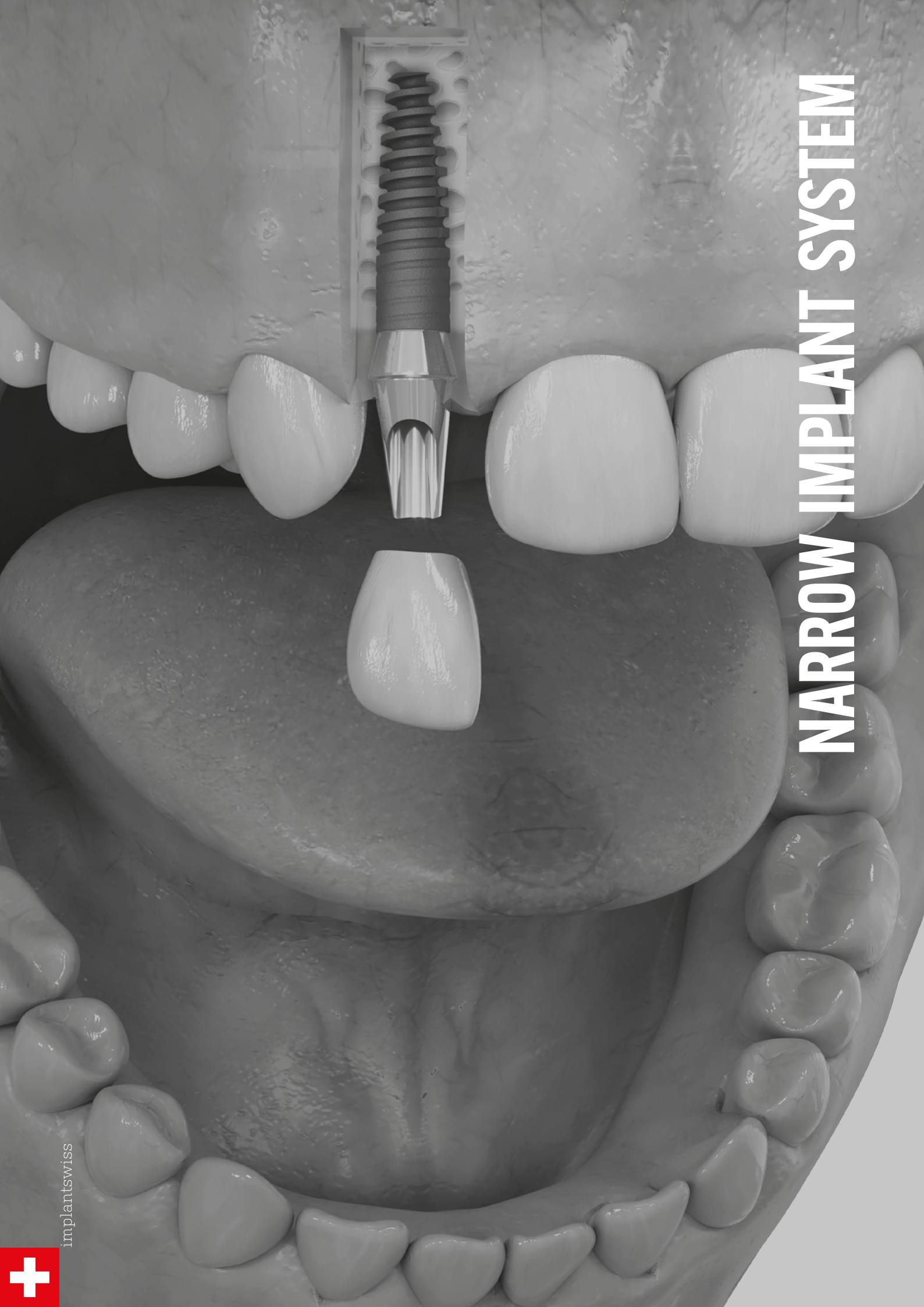

NARROW IMPLANT SYSTEM

implantswiss



**Bone Level Narrow
Hybrid Fixture**

ø2.9

10mm	S-BNFHR2910
12mm	S-BNFHR2912
14mm	S-BNFHR2914

**Bone Level Narrow
Temporary Abutment**

ø3.5

14.85mm	S-BNTA35
---------	----------

**Bone Level Narrow
Couple Abutment**

ø3.5

2.0mm	S-BNCA35207
3.0mm	S-BNCA35307
4.0mm	S-BNCA35407
6.0mm	S-BNCA35607

**Bone Level Narrow
Abutment Screw**

ø1.90

6.5mm	S-BNAS
-------	--------

**Bone Level Narrow
Closing Screw**

ø2.60

4.5mm	S-BNCLS
-------	---------

**Bone Level Narrow
Angled Abutment**

ø3.5

15°x2(Cuff)	S-BNAA35152
15°x3(Cuff)	S-BNAA35153

**Bone Level Narrow
Impression Coping Pick Up**

ø3.5

S- BNICP37

**Bone Level Narrow
Healing Abutment**

ø3.5

3mm	S-BNHA3530
6mm	S-BNHA3560

**Bone Level Narrow
Impression Coping Transfer**

ø3.5

S- BNICT37

**Bone Level Narrow
Ti-Base Abutment**

ø3.5

0.5mm	S-BNTiA05
1mm	S-BNTiA10
2mm	S-BNTiA20

**Bone Level Narrow
Lab Analog**

ø3.5

S-BNFLA






SURGICAL KIT




implantswiss




NARROW DRIVER

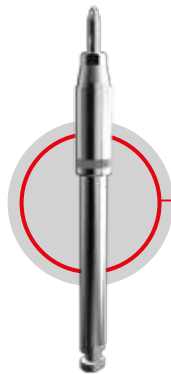

S-BNFMR

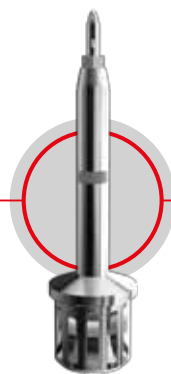

S-BNFDR

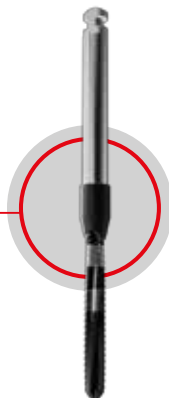

Ø2.2 LINDEMANN DRILL

S-LD2020

implantswiss.com

NOVODENT SA
Y-Parc Technopole
Avenue des Sciences 11
CH-1400, Yverdon-les-Bains
Switzerland

07.2023 / rev1

